

### L/S® Sanitary Pump Tubing Assemblies

#### Get quicker connection of tubing

##### Features/Benefits

- Meets stringent demands of biotechnology, pharmaceutical, and the food, beverage, and dairy processing industries
- Smooth, pre-molded tubing ends allow quick connection to an adapter or to another length of sanitary tubing
- Connectors help maintain cleanliness and purity of your fluids
- Sanitary tubing is less thermally conductive than metallic tubing systems and not subject to galvanic action, RFI, or EMI
- All L/S platinum-cured silicone and PharMed® BPT sanitary tubing conforms to USP Class VI and FDA standards
- FREE Certificate of resin or elastomer Compliance is available upon request

##### How to Connect Tubing

- Each length of tubing features pre-molded ends. Tubing with 1/2" mini-connections uses molded-in gasket at one end to ensure quality seal
- Pre-molded ends allow quick connection to an adapter or additional length of sanitary tubing. Order adapters below, right.
- Join 1/2" mini-connection to another 1/2" mini-connection, adapter, or to your system with a push/pull clamp

##### How to Order

1. Order your desired tubing formulation and length. Choose tubing size based on flow range (see "Specifications" on pages 700 and 702).
2. Order PVDF push/pull clamp (separately below).
3. Order silicone or Viton® gaskets as needed to connect to your system or to adapters
4. Select sanitary adapters as needed




Platinum-cured silicone tubing with pre-molded 1/2" mini connections



Tubing size	Platinum-cured silicone (96410-series) 5-ft (1.5-m) length		Platinum-cured silicone (96410-series) 10-ft (3-m) length		PharMed BPT (06508-series) 5-ft (1.5-m) length		PharMed BPT (06508-series) 10-ft (3-m) length	
	Catalog number	Price	Catalog number	Price	Catalog number	Price	Catalog number	Price
<b>L/S Precision sanitary pump tubing</b>								
L/S 13	<a href="#">GY-96100-13</a>		<a href="#">GY-96101-13</a>		<a href="#">GY-96112-13</a>		<a href="#">GY-96113-13</a>	
L/S 14	<a href="#">GY-96100-14</a>		<a href="#">GY-96101-14</a>		<a href="#">GY-96112-14</a>		<a href="#">GY-96113-14</a>	
L/S 16	<a href="#">GY-96100-16</a>		<a href="#">GY-96101-16</a>		<a href="#">GY-96112-16</a>		<a href="#">GY-96113-16</a>	
L/S 25	<a href="#">GY-96100-25</a>		<a href="#">GY-96101-25</a>		<a href="#">GY-96112-25</a>		<a href="#">GY-96113-25</a>	
L/S 17	<a href="#">GY-96100-17</a>		<a href="#">GY-96101-17</a>		<a href="#">GY-96112-17</a>		<a href="#">GY-96113-17</a>	
L/S 18	<a href="#">GY-96100-18</a>		<a href="#">GY-96101-18</a>		<a href="#">GY-96112-18</a>		<a href="#">GY-96113-18</a>	
<b>L/S High-Performance Precision sanitary pump tubing</b>								
L/S 15	<a href="#">GY-96100-15</a>		<a href="#">GY-96101-15</a>		<a href="#">GY-96112-15</a>		<a href="#">GY-96113-15</a>	
L/S 24	<a href="#">GY-96100-24</a>		<a href="#">GY-96101-24</a>		<a href="#">GY-96112-24</a>		<a href="#">GY-96113-24</a>	
L/S 35	---	---	---	---	<a href="#">GY-96112-35</a>		<a href="#">GY-96113-35</a>	
L/S 36	---	---	---	---	<a href="#">GY-96112-36</a>		<a href="#">GY-96113-36</a>	

### Accessories for L/S Sanitary Pump Tubing

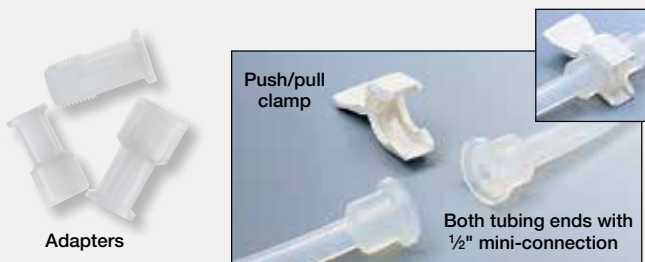

**GY-31201-88 Push/pull clamp**, PVDF. For quick joining of two 1/2" mini connections

**Gaskets.** Order gaskets to join 1/2" mini connection without molded-in gasket to another 1/2" mini connection without molded-in gasket or to an adapter (sold separately at right).

**GY-30548-00 Silicone gasket** for use with sanitary silicone (platinum-cured) tubing. Pack of 10

**GY-30548-20 Viton® gasket** for use with sanitary PharMed® BPT tubing. Pack of 10

#### Adapters

Adapter connections	Polypropylene		PVDF	
	Cat. no.	Price	Cat. no.	Price
1/2" mini to 1/8" NPT(M)	<a href="#">GY-31200-00</a>		<a href="#">GY-31201-00</a>	
1/2" mini to 1/8" NPT(F)	<a href="#">GY-31200-10</a>		---	---
1/2" mini to 1/4" NPT(M)	<a href="#">GY-31200-01</a>		<a href="#">GY-31201-01</a>	
1/2" mini to 1/4" NPT(F)	<a href="#">GY-31200-11</a>		<a href="#">GY-31201-11</a>	
1/2" mini to 1/4" hose barb	<a href="#">GY-31805-25</a>		<a href="#">GY-31808-25</a>	
1/2" mini to 3/8" NPT(M)	<a href="#">GY-31200-02</a>		<a href="#">GY-31201-02</a>	
1/2" mini to 3/8" NPT(F)	<a href="#">GY-31200-12</a>		<a href="#">GY-31201-12</a>	
1/2" mini to 3/8" hose barb	<a href="#">GY-31805-26</a>		<a href="#">GY-31808-26</a>	
1/2" mini to 1/2" NPT(M)	<a href="#">GY-31200-03</a>		<a href="#">GY-31201-03</a>	
1/2" mini to 1/2" NPT(F)	<a href="#">GY-31200-13</a>		<a href="#">GY-31201-13</a>	
1/2" mini to 1/2" hose barb	<a href="#">GY-31805-27</a>		---	---
1/2" mini to 1" maxi	<a href="#">GY-31200-30</a>		---	---
1/2" mini to 1" ladsch	<a href="#">GY-31805-05</a>		<a href="#">GY-31201-40</a>	
1/2" mini to female luer lock	<a href="#">GY-31200-50</a>		<a href="#">GY-31201-50</a>	